

AD-R169 600


A COMPUTERIZED SYSTEM FOR MEASURING DETECTION  
SENSITIVITY OVER THE VISUAL FIELD(U) ARMY RESEARCH INST  
OF ENVIRONMENTAL MEDICINE NATICK MA J L KOBRICK ET AL.  
JUN 86 USARIEM-A-42/86


1/1

UNCLASSIFIED

F/G 5/8

NL

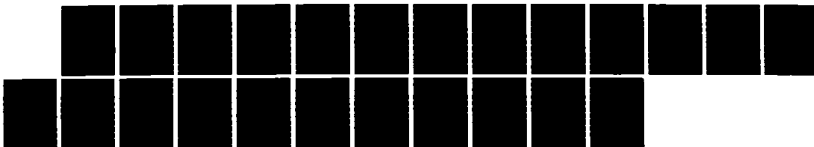

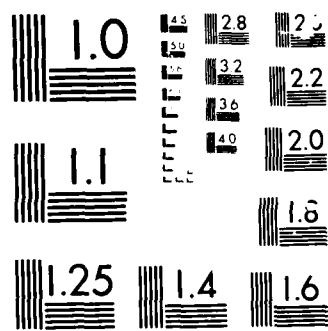

AD-A169 600

DTIC FILE COPY

UNCLASSIFIED

SECURITY CLASSIFICATION OF THIS PAGE (When Data Entered)

1

REPORT DOCUMENTATION PAGE		READ INSTRUCTIONS BEFORE COMPLETING FORM
1. REPORT NUMBER M42/86	2. GOVT ACCESSION NO.	3. RECIPIENT'S CATALOG NUMBER
4. TITLE (and Subtitle) A Computerized System for Measuring Detection Sensitivity Over the Visual Field		5. TYPE OF REPORT & PERIOD COVERED Manuscript
		6. PERFORMING ORG. REPORT NUMBER
7. AUTHOR(s) John L. Kobrick, M.S., Ph.D., Adrien R. Lussier, B.S., M.S., Stephen Mullen, B.S. and Calvin Witt, B.S., M.S.		8. CONTRACT OR GRANT NUMBER(s)
9. PERFORMING ORGANIZATION NAME AND ADDRESS US Army Research Institute of Environmental Medicine Natick, MA 01760-5007		10. PROGRAM ELEMENT, PROJECT, TASK AREA & WORK UNIT NUMBERS 3M161102BS10 64882101011
11. CONTROLLING OFFICE NAME AND ADDRESS Same as 9.		12. REPORT DATE June 1986
		13. NUMBER OF PAGES 21
14. MONITORING AGENCY NAME & ADDRESS (if different from Controlling Office)		15. SECURITY CLASS. (of this report) Unclassified
		15a. DECLASSIFICATION/DOWNGRADING SCHEDULE
16. DISTRIBUTION STATEMENT (of this Report) Approved for public release; distribution is unlimited.		
17. DISTRIBUTION STATEMENT (of the abstract entered in Block 20, if different from Report) DTIC JUL 07 1986		
18. SUPPLEMENTARY NOTES A		
19. KEY WORDS (Continue on reverse side if necessary and identify by block number) Peripheral vision, Visual field, Perimetry, Visual perception		
20. ABSTRACT (Continue on reverse side if necessary and identify by block number) A device is described for comprehensively measuring the capability of human operators to detect and react to visual signals occurring at a variety of locations throughout the functional visual field. The system is completely automated and computerized, and provides both documentation files and descriptive graphics of the operator's performance immediately upon completion of testing. Sensitivity of detection can be measured for three stimulus colors (red, yellow, green), in a testing configuration which mimics commonplace surveillance of the ordinary viewing field. An abbreviated listing of the main operating program software is provided.		

DD FORM 1 JAN 73 1473

EDITION OF 1 NOV 65 IS OBSOLETE

UNCLASSIFIED

SECURITY CLASSIFICATION OF THIS PAGE (When Data Entered)

Title:

A Computerized System for Measuring Detection Sensitivity  
Over the Visual Field

Authors:

John L. Kobrick, M.S., Ph.D.\*, Adrien R. Lussier, B.S., M.S.,  
Stephen Mullen, B.S. and Calvin Witt, B.S., M.S.

Laboratory of Origin:

US Army Research Institute of Environmental Medicine  
Natick, Massachusetts 01760

Telephone:

(617) 651-4885

Running Head:

Visual Field Tester

\*Dr. John L. Kobrick is a Research Psychologist assigned to the Health &  
Performance Division, US Army Research Institute of Environmental  
Medicine

A-1

X



Abstract

A device is described for comprehensively measuring the capability of human operators to detect and react to visual signals occurring at a variety of locations throughout the functional visual field. The system is completely automated and computerized, and provides both documentation files and descriptive graphics of the operator's performance immediately upon completion of testing. Sensitivity of detection can be measured for three stimulus colors (red, yellow, green), in a testing configuration which mimics commonplace surveillance of the ordinary viewing field. An abbreviated listing of the main operating program software is provided. *K. J. [unclear]*

Index Terms:

Peripheral vision

Visual field

Perimetry

Visual perception

One important aspect of the ongoing visual process is the detection of events both in the central and the peripheral field areas during normal surveillance of the overall visual field. Although the central line of sight involving the retinal fovea is clearly the locus of greatest visual resolution under photopic viewing conditions, the capability for detection of visual events is also present to various lesser degrees in virtually all areas of the visual field. It has been shown that the zone of best overall detection sensitivity is distributed in an ovaloid configuration about the horizontal axis of the visual field, with maximum sensitivity occurring along that axis (1). Sensitivity diminishes progressively for more peripheral locations, with the greatest losses occurring along the superior and inferior vertical meridians, probably due in part to obstructions of the visual field by the nose, chin and forehead (2,3,4,5,6,7,8).

A variety of conditions can act to degrade this basic configuration of detection capability; e.g., pathology, such as glaucoma and retinitis pigmentosa; environmental extremes, such as hypoxia; atmospheric pollution, such as carbon monoxide; and others. For these and numerous other practical reasons, it is important to be able to document accurately the zones of sensitivity for stimulus detection throughout the entire visual field, and a variety of perimeters have been developed for this purpose (e.g., Goldmann, Harms, Aimarck, Bausch & Lomb, Topcon).

Clinical perimeters provide accurate and reproducible clinical data on the state of the visual field when properly used. However, the visual task on which they are based bears little or no resemblance to commonplace viewing

in the real world. Furthermore, virtually all perimeters intrinsically, and in fact deliberately, provide cues for prior location of the test stimulus. A good example of this is the style of instrument based on use of a rotatable hemispheric surface strip upon which the stimulus is projected as a light patch. Here the subject is literally shown in advance the line of sight where the next stimulus is certain to appear. Other types of perimeters by design deliberately present a locator annulus just prior to onset of the stimulus. Thus, for clinical measurement purposes, instruments of these types alter the normal course of ordinary viewing.

The clinical measurement procedure itself is an even more general source of location cueing, which further distorts the subject's "natural" response by compartmentalizing the performance into a pre-established series of discrete observations in arbitrarily chosen segments of the visual field. Consequently, although perimeters can provide an accurate picture of the absolute limit contours of the sensate visual field, they fail to mimic the overall continuous nature of commonplace viewing, and are unable by design to assess detection sensitivity throughout the visual field area in an ongoing functional manner resembling normal vision.

The instrument described here is an automated programmable system for continuously assessing relative stimulus detection sensitivity over a variety of locations distributed throughout the visual field. This device avoids the problems of measurement bias inherent in perimeters by presenting the stimulus for detection in an unpredictable pattern of temporal occurrence and location. The subject's detection efficiency is judged on the basis of his

response time for signal detection, a measure of performance which after minimal practice quickly becomes asymptotic and remains relatively unchanged thereafter.

#### MATERIALS AND METHODS

The stimulus configuration consists of 32 units of 3 6VDC light-emitting diodes (LED's) each, distributed about the visual field of the subject. Each unit consists of one 1/8 inch diameter circular red, yellow and green LED, respectively, spaced 1/4 inch apart and mounted on a small metal template. By means of the control programming given below, this system provides the option of testing red, green or yellow sensitivity either as separate administrations or interleaved in any desired pattern.

The LED units are mounted on a hemispherically-shaped ribbed structure, which is shown schematically in Figure 1.

-----  
Figure 1 about here  
-----

This hemisphere is constructed of 3/32 inch x 3/4 inch curved aluminum strips, each joined at one end to an 8-inch diameter circular center-plate at equal intervals around its edge, and connected sequentially at the other end by 19-inch spacing rods. This arrangement results in a structure composed of eight curved ribs spaced at 45-degree intervals to form a 48-inch diameter hemisphere. The interior surfaces of the ribs are faced with the black pile membrane of hook-and-pile fastener material (Velcro), and a tab of the accompanying hooked membrane is cemented to the back of each LED template, so



that the stimulus lights can be conveniently adhered anywhere along the ribs of the hemisphere. In the research in which the apparatus has currently been employed, the LED units were positioned on each of the eight ribs at  $12^{\circ}$ ,  $38^{\circ}$ ,  $64^{\circ}$  and  $90^{\circ}$  angular displacement from the center of the hemisphere. However, the Velcro feature allows the positioning of the LED stimuli in any desired arrangement. It should also be noted that the common radius of the hemispheric rib structure provides that the stimulus lights are all at the same distance from the viewing position, and thus subtend the same visual angle.

All metal surfaces of the rib structure are painted flat black, and the back side of the entire assembly is enclosed in a matte black fabric cover. The center plate of the hemisphere is mounted by a sliding clamp to a vertical pipe to allow adjustment of viewing height to the subject's seated position.

The subject is seated in a contoured seat with an adjustable headrest. (A salvaged Volkswagen front bucket seat with an adjustable reclining control was found to be cheap and effective.) The seat is positioned so that the subject's line of sight intersects the center of the hemisphere, and the edge of the hemisphere abuts the peripheral limit ( $90^{\circ}$ ) of the visual field.

The entire task is generated and administered by a Hewlett-Packard HP-07XM computer and an associated H-P Model 82901M dual floppy-disk drive. The following HP ROM's are also required: I/O; Plotter; Advanced Programming; Matrix; and, 256 Kbytes Memory. This system is configured in turn with two H-P Model 3495A scanners, each with three 20-channel low thermal relays and

one 10-channel relay actuator. These latter units generate the necessary switch closures for LED activation. Additional devices used in the present assembly include an H-P Model 59313A A/D converter, and an H-P Model 9871A character impact printer. A schematic diagram of this equipment system is shown in Figure 2.

-----  
Figure 2 about here  
-----

#### Operating Procedure and Software

The subject's task is to orient his line of sight toward the center of the display, and to monitor the entire visual area for activations of any of the LED's. He is trained to respond to any detection by pressing a push-button switch held in the preferred hand. The system measures and stores the response time interval (RT) in milliseconds following each LED activation to the point of switch closure. Thus, RT for signal detection is the basic measure employed by this system as an index of relative visual field sensitivity. As presently programmed, the system presents all lights of the specified class (red, green, yellow or their combinations) in randomized order of occurrence, with the provision that no light is repeated until all of the specified class have occurred. The lights are presented at randomized intervals ranging in value between 5 and 20 seconds following switch closure of the previous detection. Signals which are not detected within 5 seconds are considered to be missed, and the system then arbitrarily assigns them that RT value, which is a virtually infinite magnitude in the realm of

response time or reaction time behavior, and, therefore, a missed signal. In practice, the RT has proven to be a good measure to represent relative visual detection capability, because the task is acquired quickly, remains asymptotic once training is complete, and appears to be sensitive to exposure to environmental and other stressors.

As presently configured, all system software is programmed in Hewlett-Packard BASIC. A complete copy of the software can be obtained by contacting the authors.

After boot-up and initialization, the program prompts for insertion of system and data disks and checks for correct insertion. It then prompts for calibration of the A/D converter and for initiation of output devices (printer or other accessories). (A/D calibration is recommended at least once daily, or every four hours when the system is used continuously.) Prompts are then given for identification of the test conditions, stimulus characteristics, and historical documentation of the test run and research study involved. The system then requests identification of the test subject (name, social security number, optional information) and identification of the technician. All instructions thereafter are user-friendly, in that the technician is referred to by name. The program includes numerous reminders and helpful prompts to assist in operation of the equipment. Key-stroke sub-routines are also provided for performing system checks of the stimulus lights, in which all lights flash individually in sequence.

In its present form, the program provides a READY prompt, and when given a RUN command, conducts the entire testing procedure according to the

stimulus and test parameters specified at the beginning of the procedure. Instructions are provided to the subject automatically through a voice synthesizer at the appropriate points in conduct of the test. Upon completion of the testing procedure, the system stores the RT data and provides a printed listing of the relevant test information, and the RT values listed against the stimulus characteristics and visual field locations. Finally, the system prints two graphic representations of the averaged RT's, one for RT versus the degree of peripheral stimulus location, and the other for RT versus radial position of stimuli in the visual field.

#### Abbreviated Software Program

An abbreviated version of the main operating program of this system is as follows.

```

10 | PERIPHERAL VISION TESTING SYSTEM
20 | USARIEM, HEALTH & PERFORMANCE DIVISION
30 |   by SFC AR LUSSIER -- VERSION 3.03 -- DECEMBER 1983
40 |
50 |
3010 | ***** LIGHTS/COLOR SECTION *****
3020 | *****
3021 |
3022 | This section prepares the light sequences for presentation
3023 | ----- User selection of light stimuli -----
3024 |
3030 OFF CURSOR @ ALPHA SCREEN @ FOR A=1 TO 10 @ DISP @ NEXT A @ OFF CURSOR @
ALPHA SCREEN @ DISP "LIGHT Information"
```

```

3040 FOR A=1 TO 3 @ ON A GOSUB Red ,Green ,Yellow @ DISP USING "4X,D,X,K" ;
A,"= "&Temp1$ @ NEXT A

3050 DISP TAB (5);"R = Random selection of color(s)" @ DISP TAB (3);"Rxx =
Random order of any two specified colors"

3060 L1=1 @ Prompt$="Number of colors to be processed" @ GOSUB Question

3070 IF LEN (Ans$)#1 OR Ans$[1,1]<"0" OR Ans$[1,1]>"3" THEN DISP "Please
enter a number within the range of 1 to 3" @ OFF CURSOR @ ALPHA I2 @ GOTO
3060

3080 IF Stdtype$="CD" THEN Color=1 ELSE Color=VAL (Ans$)

3090 IF Stdtype$="PV" THEN 3120

3100 OFF CURSOR @ ALPHA I2+2 @ DISP @ OFF CURSOR @ ALPHA I2 @ L1=1 @
Prompt$="Primary color. If you have no preference type in 'R'" @ GOSUB
Question

3110 GOTO 3130

3120 OFF CURSOR @ ALPHA I2+2 @ DISP @ OFF CURSOR @ ALPHA I2 @ L1=4 @
Prompt$="Color order. If you have no preference type in 'R'" @ GOSUB
Question

3130 IF LEN (Ans$)>4 THEN OFF CURSOR @ ALPHA I2 @ GOTO 3090

3140 Order$=UPC$ (Ans$) @ IF POS (Order$,"R")=0 AND LEN (Order$)>Color THEN
OFF CURSOR @ ALPHA I2 @ GOTO 3090

3150 IF POS (Order$,"R")#0 OR Color=1 AND POS (Order$,"R")=0 AND LEN
(Order$)=1 THEN Ans$="N" @ GOTO 3180

3160 L1=1 @ Prompt$="Repeat the sequence of lights? Y/N" @ OFF CURSOR @
ALPHA I2 @ GOSUB Question

```

```

3170 IF LEN (Ans$)=1 AND (POS (UPC$ (Ans$),"Y")=0 OR POS (UPC$ (Ans$),"N")=0)
THEN 3180 ELSE OFF CURSOR @ ALPHA I2 @ GOTO 3160
3180 Rep$=UPC$ (Ans$) @ OFF CURSOR @ ALPHA I2 @ DISP @ DISP HGL$ ("Please be
patient while I try to generate some random numbers!") @ OFF CURSOR
3190 IF Stdtype$="CD" THEN GOTO 5000
3191 !
4000 ! ***** RANDOM PROCESS VISUAL PERCEPTION *****
4010 ! =====
4011 !
4020 RANDOMIZE DATE *TIME @ Lights$,Delay$="" @ IF POS (Rep$,"Y")#0 OR POS
(Order$,"R")=0 THEN Color1=1 ELSE Color1=Color
4030 Temp2$,Temp3$="" @ FOR A=1 TO 32*Color1
4040 Colors$=CHR$ (INT (RND *32*Color1)+1) @ IF POS (Temp2$,Colors$)#0 THEN
4040 ELSE Temp2$=Temp2$&Colors$ @ Temp3$=Temp3$&CHR$ (INT (RND *10)+1)
4050 NEXT A @ IF POS (Rep$,"Y")#0 THEN Lights$=RPT$ (Temp2$,Color) @
Delay$=RPT$ (Temp3$,Color) @ GOTO 4080
4060 IF POS (Order$,"R")=0 .AND LEN (Lights$&Temp2$)/32=Color OR POS
(Order$,"R")#0 THEN 4070 ELSE Lights$=Lights$&Temp2$ @ Delay$=Delay$&Temp3$ @
GOTO 4030
4070 Lights$=Lights$&Temp2$ @ Delay$=Delay$&Temp3$
4080 FOR A=1 TO LEN (Lights$)
4090 IF NUM (Lights$[A,A])>32 THEN Lights$[A,A]=CHR$ (NUM (Lights$[A,A])-32)
@ GOTO 4090
4100 NEXT A @ Colors$="" @ MAT Results=ZER @ IF POS (Order$,"R")#0 THEN 4120

```

```

4110 FOR A=1 TO Color @ Colors$=Colors$&RPT$ (Order$[A,A],32) @ NEXT A @ GOTO
6010
4120 Temp2$,Temp3$="" @ IF Color=1 THEN Order$=VAL$ (INT (RND *3)+1) @ GOTO
4110
4130 IF LEN (Order$)#1 THEN 4170
4140 FOR A=1 TO Color
4150 Temp3$=VAL$ (INT (RND *3)+1) @ IF POS (Temp2$,Temp3$)=0 THEN
Temp2$=Temp2$&Temp3$ ELSE GOTO 4150
4160 NEXT A @ Order$="R"&Temp2$
4170 Order$=ROTATE$ (Order$,LEN (Order$)-POS (Order$,"R")+1)
4180 A=6 @ FOR B=2 TO LEN (Order$) @ A=A-VAL (Order$[B,B]) @ NEXT B @ IF A=0
THEN 4200
4190 FOR B=0 TO 32 @ Results(A,B)=1 @ NEXT B
4200 FOR A=1 TO LEN (Lights$)
4210 B=INT (RND *3)+1 @ IF Results(B,NUM (Lights$[A]))=1 THEN 4210 ELSE
Colors$[A,A]=VAL$ (B) @ Results(B,NUM (Lights$[A]))=1
4220 NEXT A @ GOTO 6010
4221 !
5000 ! =====* RANDOM PROCESS COLOR DISCRIMINATION *=====
5010 ! =====
5011 !
5020 IF Order$="R" THEN Order$=VAL$ (INT (RND *3)+1)
5030 Temp1$,Temp2$="" @ FOR A=1 TO 32
5040 Temp1$=CHR$ (INT (RND *32)+1) @ IF POS (Temp2$,Temp1$) THEN 5040 ELSE
Temp2$=Temp2$&Temp1$

```

```

5050 NEXT A @ Lights$=""
5060 B=1 @ FOR A=1 TO 32 @ Temp=INT (RND *3)
5070 Lights$[B]=RPTS ("|",Temp)&Temp2$[A,A] @ Colors$[B]=RPTS
("0",Temp)&Order$
5080 B=B+Temp+1 @ NEXT A @ B=B-1
5090 Delay$="" @ FOR A=1 TO B @ Delay$=Delay$&CHR$ (INT (RND *10)+1) @ NEXT A
5100 C=3*LEN (Delay$) @ Lights$=Lights$[1,B]&RPTS ("|",C-B) @
Colors$=Colors$[1,B]&RPTS ("0",C-B)
5110 Temp1$,Temp2$="" @ C=B @ FOR A=1 TO LEN (Delay$)
5120 Temp=INT (RND *3) @ IF Temp=0 AND NUM (Lights$[A,A])=0 THEN 5120
5130 IF Temp#0 THEN 5140 ELSE Lights$[C+A,C+A],Lights$[C*2+A,C*2+A]="| " @
Colors$[C+A,C+A],Colors$[C*2+A,C*2+A]="0 " @ GOTO 5170
5140 FOR B=1 TO Temp @ Temp1$=CHR$ (INT (RND *32)+1)
5150 Temp2$=VAL$ (INT (RND *3)+1) @ IF POS (Colors$[1,C],Temp2$)#0 OR
Colors$[C+A,C+A]=Temp2$ THEN 5150
5160 Lights$[C*B+A,C*B+A]=Temp1$ @ Colors$[C*B+A,C*B+A]=Temp2$ @ NEXT B
5170 NEXT A
5180
"
10000 ! =====
10001 !
10010 ! DATA COLLECTION
10020 ! SECTION #1 ** SET-UP
10021 !
10030 ! =====

```



```

10040 REDIM Results(0,LEN (Delay$))@  FOR RUNS=SUBJ+1 TO SUBJ+Subnum ! OUTER
SUBJECT LOOP

10050 PAGESIZE 16 @ ON KEY# 1,"START" GOTO 10110 @ OFF KEY# 2 @ OFF KEY# 8

10060 CLEAR @ MAT Results=ZER @ DISP "Subject #"&VAL$ (RUNS-SUBJ)&" is
"&NAME$(RUNS) @ DISP

10070 L1=32 @ Prompt$="Any special comment to be made?  (If none type in
NONE)" @ GOSUB Question

10080 IF UPC$ (Ans$)="NONE" THEN CMT$(RUNS)=" " ELSE CMT$(RUNS)=Ans$

10090 CLEAR @ KEY LABEL @ DISP "Press "&HGL$ ("KEY #1")&" to start." @ DISP @
DISP "Subject:  "&HGL$ (NAME$(RUNS)) @ OFF CURSOR

10100 GOTO 10100

10110 ALPHA 1 @ ALPHA @ ON KEY# 8,"ABORT" GOTO Abort @ ON KEY# 1,"RESTART"
GOTO Restart

10120 Temp1$="Subject$:  "&HGL$ (NAME$(RUNS)) @ OFF CURSOR @ ALPHA 1,INT
((80-LEN (Temp1$))/2) @ AWRT Temp1$

10130 !

10140 ALPHA 2,18 @ AWRT "Seq#  Light  Color  Delay  Reaction  Mean"

10150 RESTORE 10160 @ FOR C=8 TO 16 @ READ Temp1$,X@ ALPHA C,32+X @ AWRT
Temp1$ @ NEXT C

10160 DATA ! CONTAINS ASCII CODES FOR LIGHTS

10170 DATA ! CONTAINS ASCII CODES FOR LIGHTS

10180 DATA ! CONTAINS ASCII CODES FOR LIGHTS

10190 ALPHA 15,1 @ AWRT HGL$ (" ABORT ") @ ALPHA 16,1 @ AWRT HGL$
("RESTART") @ ALPHA 1,1 @ ON CURSOR

10191 !

```

```

11000 ! SECTION #2 ** INFORMATION DISPLAY
11010 ! -----
11011 !
11020 MAT Results=ZER @ Smean,Snum=0 @ Tmm$(RUNS)=TIMES @
Tmm$(RUNS)=Tmm$(RUNS)[1,5] @ FOR B=1 TO LEN (Delay$) @ Buf$="" @ ON EOT 7
GOTO 12050
11030 I2=4 @ FOR C=4 TO 6 @ OFF CURSOR @ ALPHA C @ DISP @ NEXT C
11040 FOR C=B-1 TO B+1 @ CFLAG 1 @ IF C<1 THEN C=B
11050 IF C>LEN (Delay$) THEN 11160
11060 Temp1$="" @ Temp1$=VAL$ (C) @ Temp1$=RPT$ (" ",4-LEN (Temp1$))&Temp1$ @
OFF CURSOR @ ALPHA I2,18 @ AWRIT Temp1$
11070 Temp1$=VAL$ (NUM (Lights$[C,C])) @ Temp1$=RPT$ (" ",5-LEN
(Temp1$))&Temp1$ @ ALPHA I2,24 @ AWRIT Temp1$
11080 X=VAL (Colors$[C,C]) @ IF X=0 THEN Temp1$="None" @ GOTO 11100
11090 ON X GOSUB Red ,Green ,Yellow
11100 ALPHA I2,31 @ AWRIT Temp1$&" "
11110 Temp1$=VAL$ (NUM (Delay$[C,C])) @ Temp1$=RPT$ (" ",5-LEN
(Temp1$))&Temp1$ @ ALPHA I2,39 @ AWRIT Temp1$
11120 IF Results(0,C)=0 THEN ALPHA I2,47 @ AWRIT RPT$ (" ",16) @ GOTO 11160
11130 Temp1$="" @ Temp1$[1,8]=VAL$ (Results(0,C)) @ ALPHA I2,47 @ AWRIT
Temp1$
11140 Temp1$="" @ IF Smean=0 OR Snum=0 THEN Temp1$="0" ELSE Temp1$[1,5]=VAL$
(Smean/Snum)
11150 ALPHA I2,58 @ AWRIT Temp1$
11160 I2=I2+1 @ IF FLAG (1) THEN RETURN

```

```

11170 IF C>LEN (Lights$) THEN 11190
11180 NEXT C
11190 IF I2=6 AND (B=1 OR B=LEN (Lights$)) THEN I2=4
11200 IF I2=6 OR I2=7 THEN I2=5
11210 CFLAG 1 @ GOSUB Onlight
11211 !
12000 ! SECTION #3 ** CHANNEL SELECTION AND DATA COLLECTION
12010 ! -----
12011 !
12020 OUTPUT 717 ;"H4BJ" @ OUTPUT 718 USING "#,K" ; OUT$(B)
12030 X=TIME @ TRANSFER 717 TO Buf$ INTR ; DELIM 2 COUNT 1590
12040 GOTO 12040

12050 X1=TIME @ OUTPUT 718 ;"C" @ IF X1-X>= 5 THEN X,Results(0,B)=5 ELSE
X,Results(0,B)=X1-X
12060 IF X#5 THEN Smean=Smean+X @ Snum=Snum+1
12070 ON TIMER# 1,1000*NUM (Delay$(B,B))+5-X GOTO 12100
12080 SFLAG 1 @ C=B @ GOSUB 11050
12090 GOTO 12090
12100 OFF TIMER# 1 @ IF Stdtype$="PV" THEN GOSUB Offlight ELSE GOSUB Multoff
12110 NEXT B

```

## REFERENCES

1. Kobrick JL. Effects of physical location of visual stimuli on intentional response time. J. Eng. Psychol. 1965; 4:1-8.
2. Kobrick JL. Effects of hypoxia on response time to peripheral visual signals. In The perception and application of flashing lights, Proceedings of an international symposium. London: Hilger, Ltd, 1971.
3. Kobrick JL. Effects of hypoxia on voluntary response time to peripheral visual stimuli during central target monitoring. Ergonomics. 1972; 15:147-156.
4. Kobrick JL. Effects of hypoxia on peripheral visual response to rapid sustained stimulation. J. Appl. Physiol. 1974; 37:75-79.
5. Kobrick JL. Effects of hypoxia on peripheral visual response to dim stimuli. Perceptual and Motor Skills. 1974; 41:467-474.
6. Kobrick JL, Appleton B. Effects of hypoxia on visual performance and retinal vascular state. J. Appl. Physiol. 1971; 31:357-362.
7. Kobrick JL, Dusek ER. Effects of hypoxia on voluntary response time to peripherally placed visual stimuli. J. Appl. Physiol. 1970; 29: 444-448.
8. Kobrick JL, Sutton WR. Device for measuring voluntary response time to peripherally placed stimuli. Percep. Mot. Skills. 1970; 35: 255-258.

Disclaimer

The views, opinions and/or findings contained in this report are those of the author(s), and should not be construed as an official Department of the Army position, policy or decision, unless so designated by other official documentation.



Human subjects participated in these studies after giving their free and informed voluntary consent. Investigators adhered to AR 70-25 and USAMRDC Regulation 70-25 on Use of Volunteers in Research.


Citations of commercial organizations and trade names in this report do not constitute an official Department of the Army endorsement or approval of the products or services of these organizations.

Figure Captions

Figure 1. Schematic diagram of the visual testing system

Figure 2. Schematic diagram of equipment system

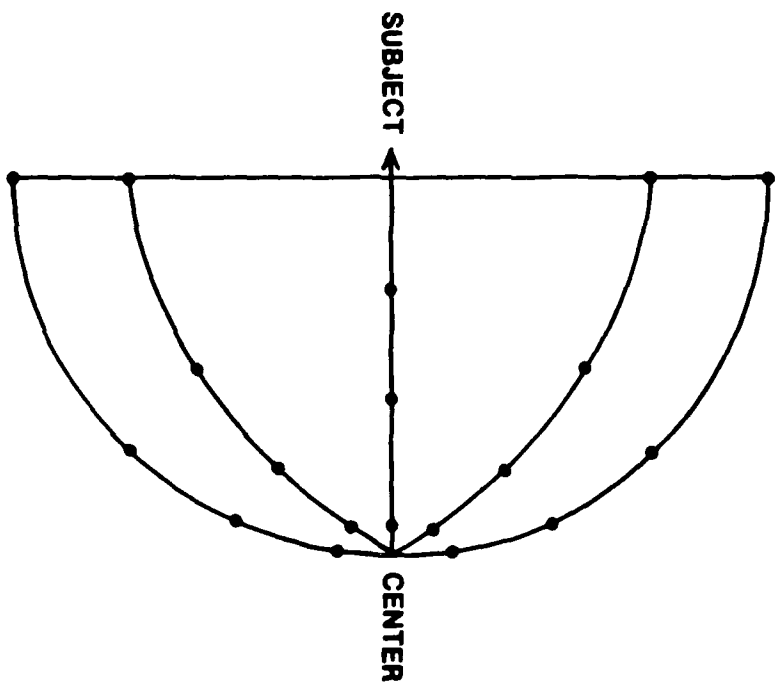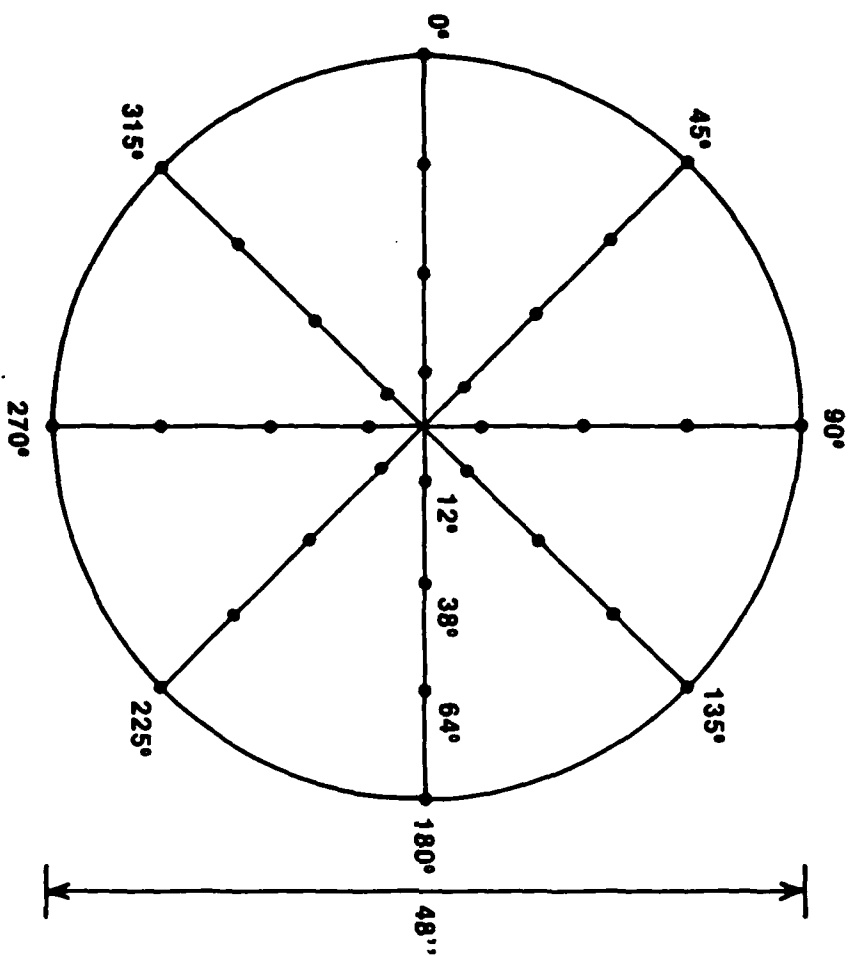

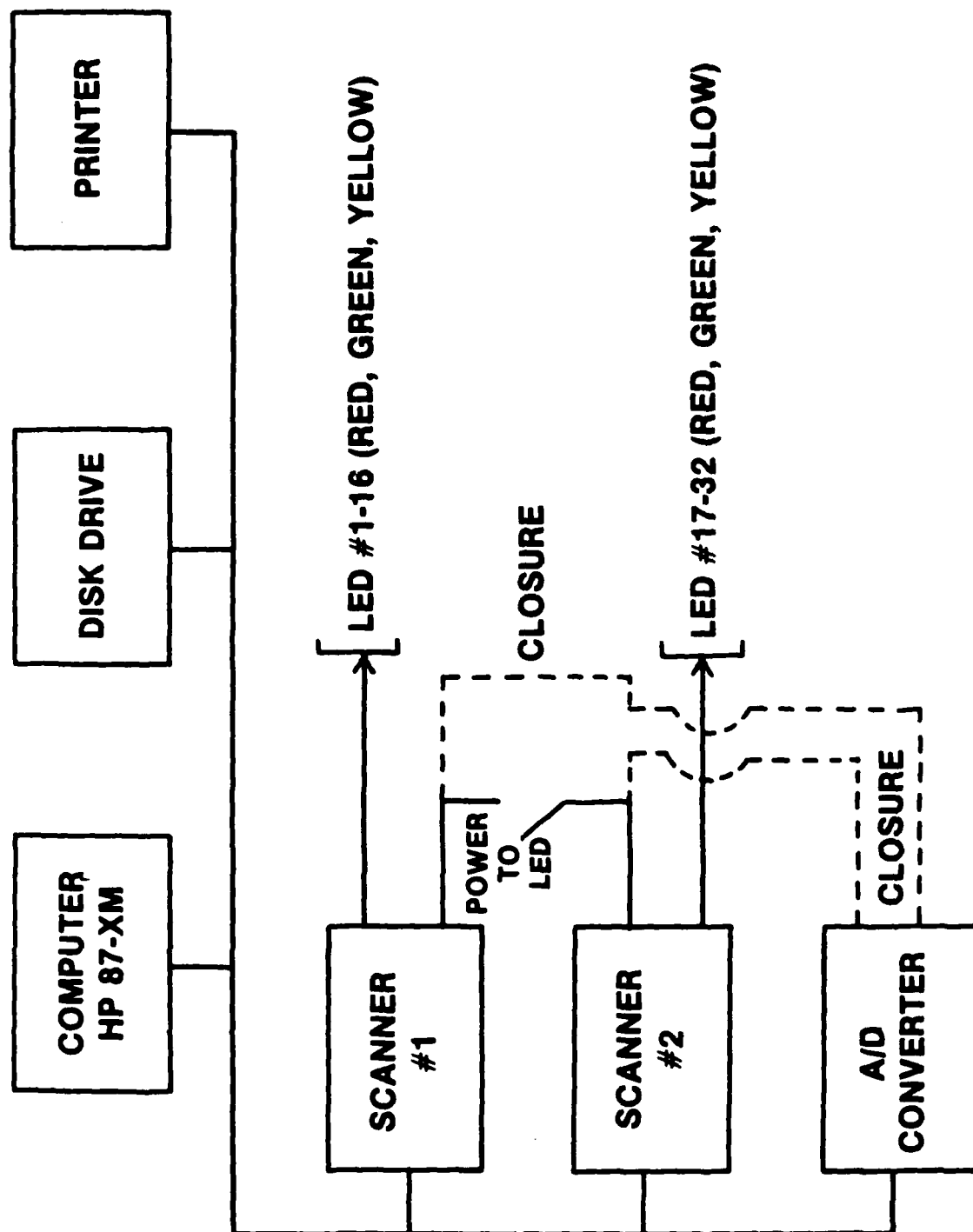



END

DT/C

8-86